










DS 360

Fully Automated Hemoglobin A1c Testing

Building on DREW Scientific's extensive experience in developing Diabetic Testing Systems. The DS360 is a fast, easy to use HbA1c Analyser.

-  Based on HPLC technology
-  A1c results in minutes
-  Offers immediate risk assessment for diabetic complications
-  Patient reports printed automatically
-  A1c Quantitation in the presence of HbS, HbC and HbF
-  Optimized to minimize interference from carbamylation, lipemia and labile A1c
-  1000 tests reorder pack
-  Calibrate A1c once/cartridge
-  Update via USB port A1c program parameters







The DS360 continues the long line of diabetic testing products from Drew Scientific. Drew Scientific has been developing diabetic systems since the mid 1980's and the DS360 is the culmination of this extensive experience.

The DS360 uses the reference method of ion-exchange HPCL to accurately determine A1c in minutes. The same technology is used to quickly identify Hemoglobin Variants.

The DS360 provides a fully automated system with on board patient demographic database management and QC capabilities. The software is an easy to use Graphical User Interface and has multiple language options available. The DS360 will be traceable to the IFCC reference method and will be NGSP Certified.

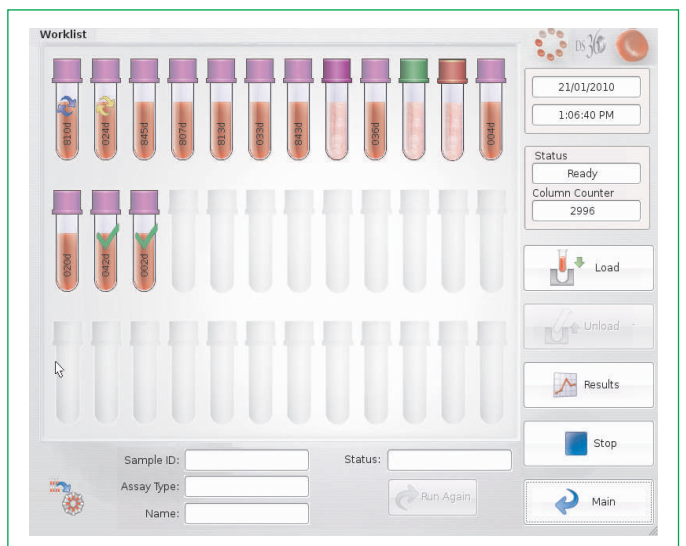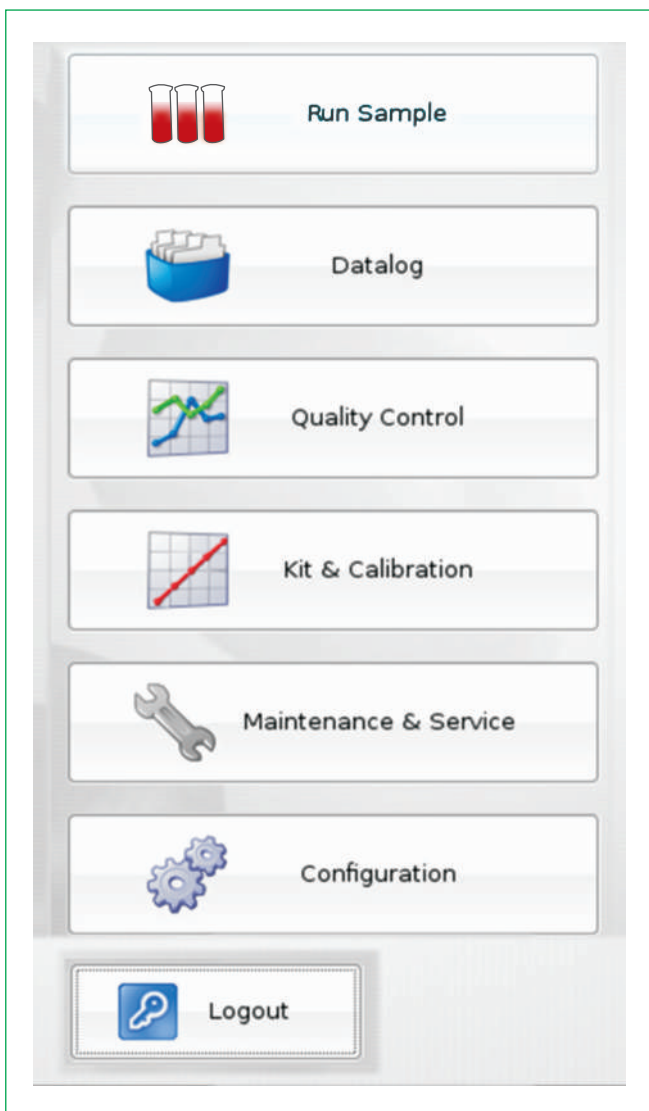

Exclusively Marketed in India by :

TRANSASIA BIO-MEDICALS LTD.

Transasia House, 8 Chandivali Studio Road, Andheri (E), Mumbai - 400 072

Tel.: (022) 40309000, Fax : (022) 2857 3030 Email : transasia@transasia.co.in Website : www.transasia.co.in



Drew Scientific Inc is an ISO registered company
Drew Scientific is a division of...



Drew Scientific Inc © 2007